




The SHOEBBOX Solution


Testing Equipment. Data Management. Professional Services. Customer Care.

SHOEBOX Audiometry is the world's first automated iPad audiometer. More than simply portable, SHOEBOX performs diagnostic hearing testing using a modified Hughson-Westlake protocol. It is optimized—and validated—for use outside of a sound booth thanks to its patent pending Response and Environment Adaptive Control Technology (REACT™). REACT ensures test accuracy by continuously monitoring patient response behavior and environmental conditions such as ambient noise levels. It alerts administrators when it detects anomalies that may impact results accuracy, enabling immediate corrective action.

Robust Functionality


- ✓ Manual, automated, and assisted testing modes
- ✓ REACT Wizard
- ✓ Portable: iPad, calibrated headphones, software, web-based data management
- ✓ Pure tone: air conduction with masking (Standard), air & bone with masking (Pro)
- ✓ Speech reception threshold (SRT)
- ✓ Speech discrimination testing (NU-6, PBK word lists)
- ✓ Automated shift detection
- ✓ Extended high frequencies (to 16kHz)
- ✓ Embedded inventories and custom questionnaires

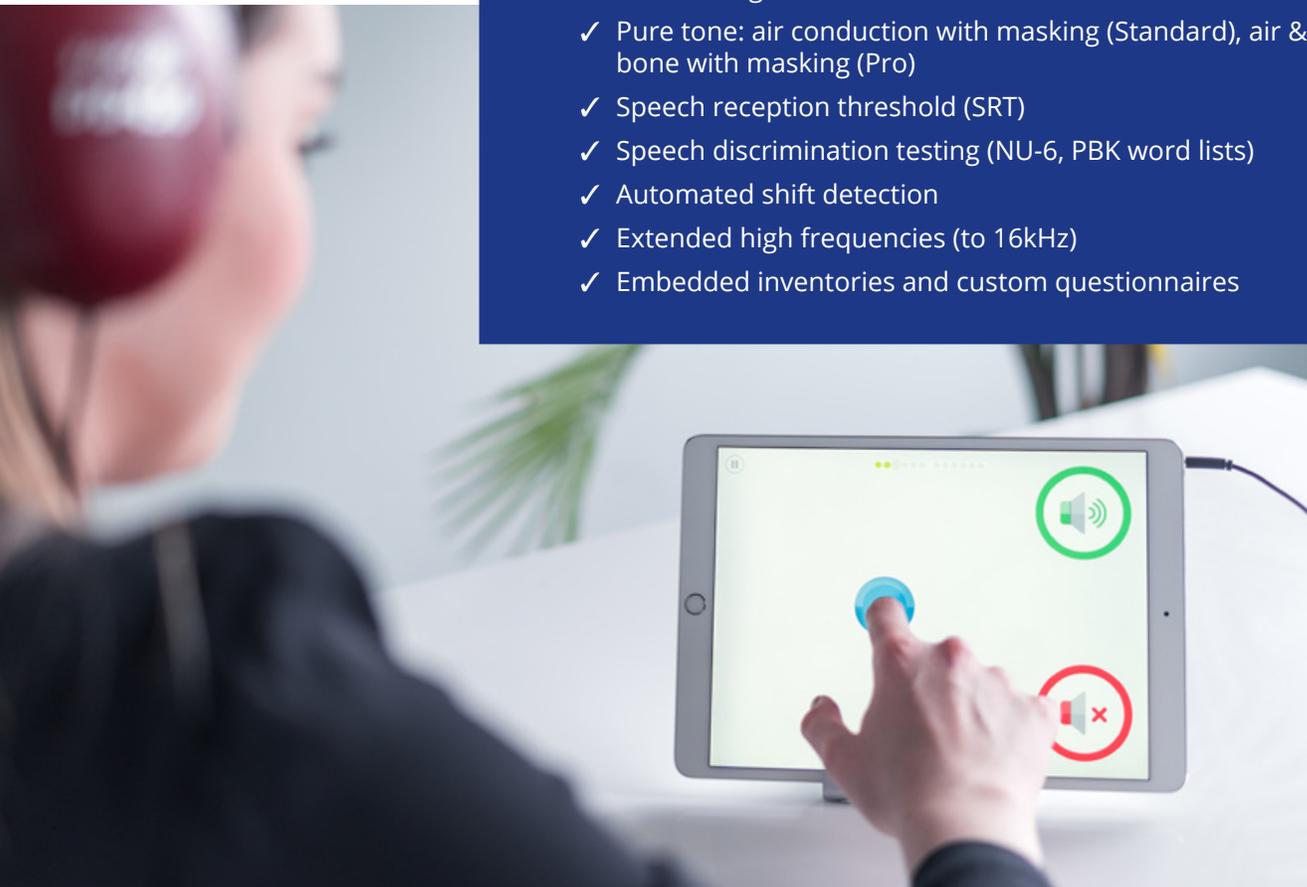

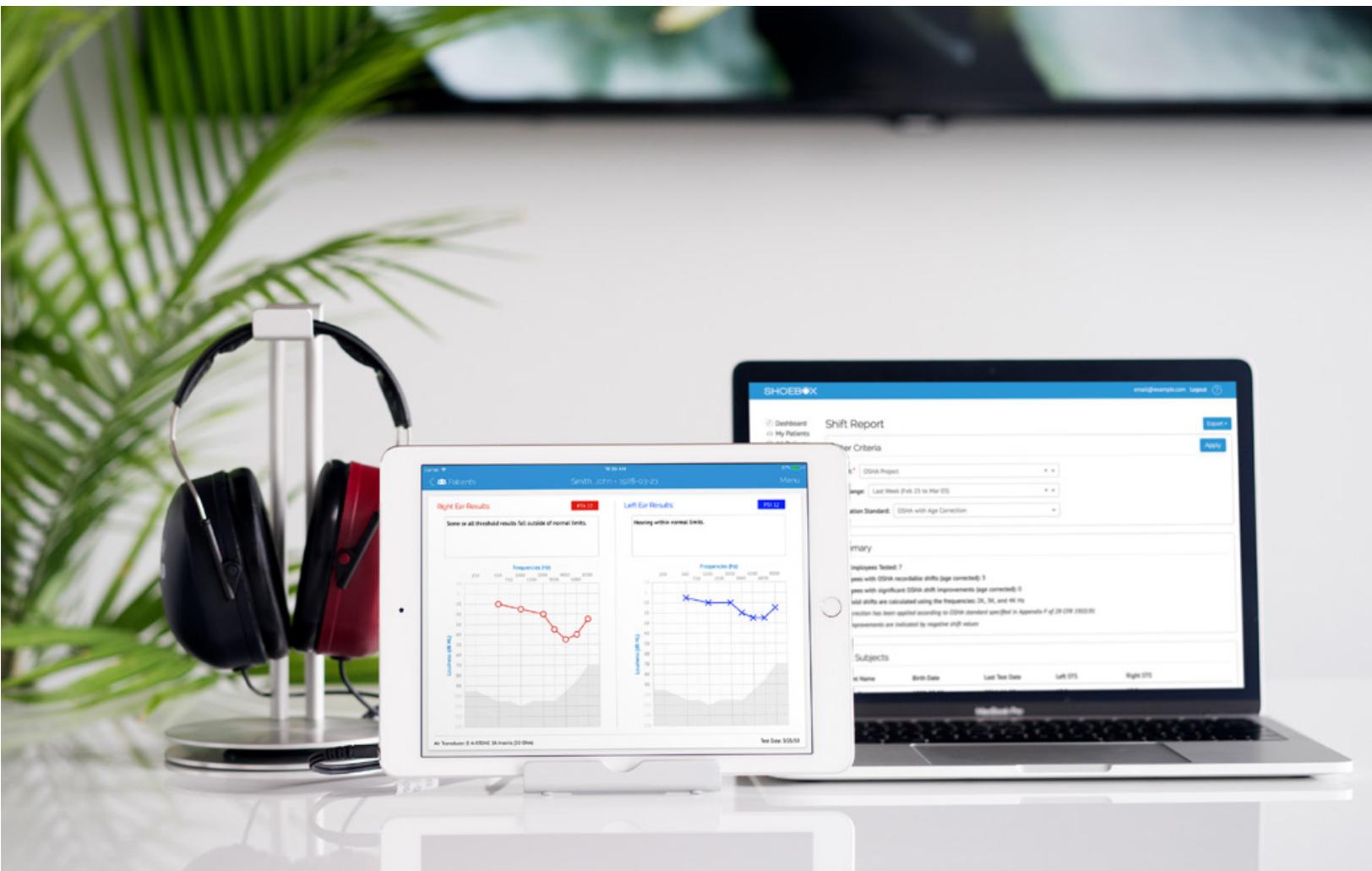

Easy-to-Use and Portable

Comprised of the testing equipment, a web-based data management system, professional services, and Customer Care, SHOEBOX is the easiest to use and most portable audiometric system available on the market today. Ideal for mobile clinicians, physicians, clinical researchers, occupational hearing conservationists, and more, SHOEBOX packs all the power of a traditional audiometer into a small, mobile form factor. The system offers a variety of flexible testing options, from self-testing to test administration by hearing health professionals and even non-audiologists.


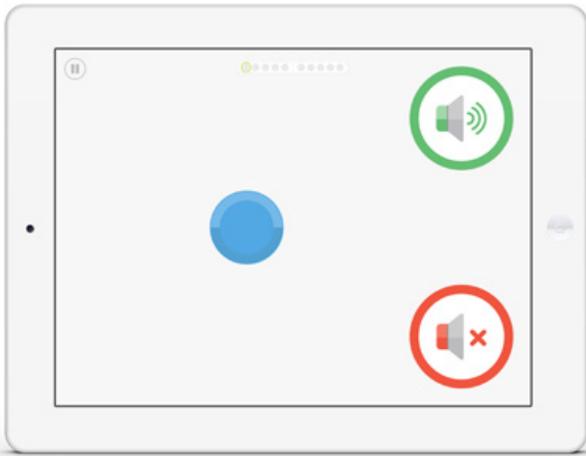


Automated Game Play

The user simply taps the “sound disc” and drags it to the red or green targets based on whether or not they hear various tones presented at different frequencies. The user can also tap the targets to answer instead of dragging.

Assisted Mode

Gives the test administrator an option for faster testing of patients who may have challenges with the automated test.

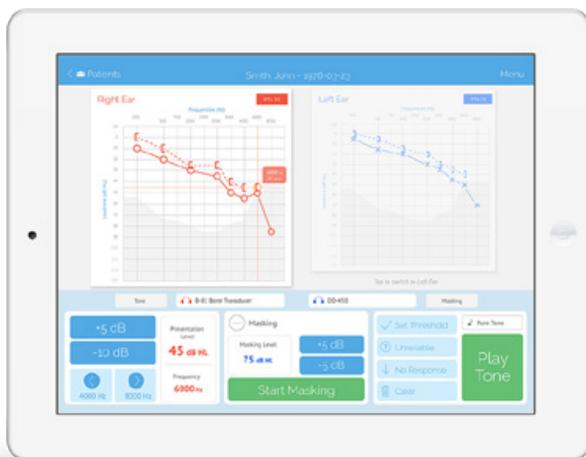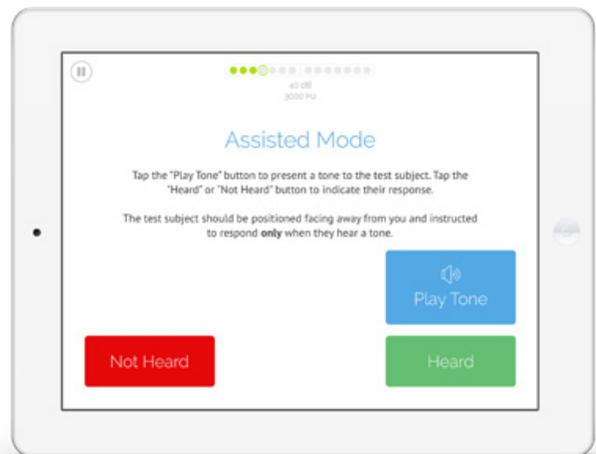


Full Manual Control

Preferred by audiologists, ENT physicians, and CAOHC certified hearing conservationists, SHOEBOX offers full manual control with advanced clinical capabilities. It can also be used to perform functional testing with select cochlear implant sound processors.


Data Management, Filtering, and Reporting Capabilities

The SHOEBBOX Data Management platform is unmatched in its ability to manage and analyze patient demographic data and test results across various projects, sites, and SHOEBBOX audiometers. It offers reporting capabilities designed to help you adhere to your specific regulatory compliance requirements. Our Data Management system is a HIPAA-compliant web-based portal used to securely archive data and audiograms, speech tests, and questionnaires captured on the iPad.


If you manage large-scale projects—or have numerous clients, employees, patients, or even clinical trial subjects—Data Management provides an environment where information can be stored, organized, and viewed in meaningful ways, and where access privileges can be granted or removed as employees shift between projects. Extensive group and individual summary reports can be generated and provided for review based on a number of threshold shift and classification grade criteria. Whether you test patients at a single clinic, clients or employees at multiple locations for hearing conservation programs, or even if you are involved in hearing testing for research or investigation purposes, our solution provides the tools you need to gather and analyze your data.

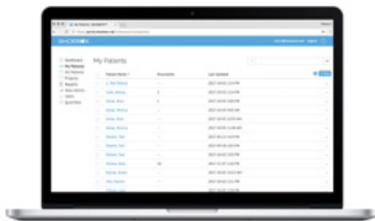

View and filter patient lists

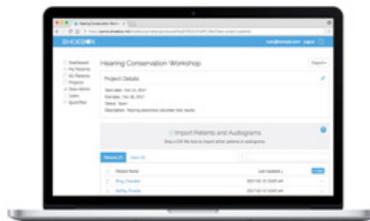

Organize data into logical groupings

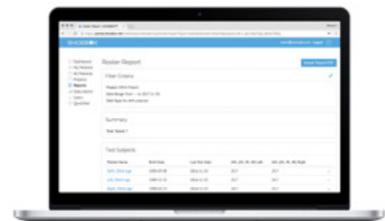


Roster Report: Summary list of individuals tested in a given time period



Shift Reports: List of individuals where shift criteria has been detected



Individual patient/employee report with consolidated data from all tests and questionnaires

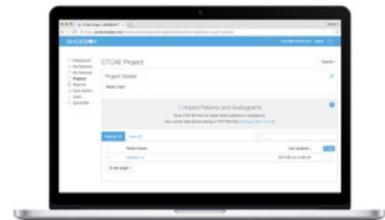

Data Import and Export functionality


Customer Care and Professional Services

Not only are we the manufacturers of SHOEBOX Audiometry, but we also service, train, and support our customers. Our goal is to ensure your 100% satisfaction and success with SHOEBOX. Our Customer Care and Professional Services include:



Training

We include 1 hour of “Getting Started” training and access to our online “Audiometry Essentials” training modules. Contact us if you require further customized training or just wish to learn more.

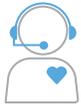

Technical Support

Whether by phone, email, web, or video chat, our support staff members are here to answer all of your questions and provide you with expert technical assistance.



Success Check-ins

Once you are up and running, our Customer Care team will continue to reach out and check in. These success check-ins provide an opportunity to ask questions, give feedback, and make requests. We value these personal interactions and believe they will help you be more successful with your SHOEBOX solution.



Device Management

We remotely manage the iPads for you. From account creation to software updates, this service is completely turnkey. When the systems arrive, just turn them on and you're ready to go! We take care of everything. We can even remotely wipe them for you if they need to be deployed to new users.

Annual Calibration Service

To maintain compliance with ANSI S3.6 standards, we offer:



- Annual notification when calibration is due
- Automatic swap out of newly calibrated transducers
- New calibration certificate for ANSI records management

We also offer customized services tailored to your needs. Some examples include:




Deployment Services

Regardless of how many systems you need and no matter where in the world you need them delivered, we have the expertise to handle all of your deployment needs.



Translations

SHOEBOX Audiometry is currently available in English, French, Spanish, and Portuguese. The patient-facing screens of the automated test are available in 27 languages. If you need a language we don't currently have, we offer custom translation services.

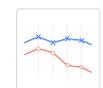

Audiological Consulting

Our in-house licensed audiology team provides support and consultation when you need it. Whether for protocol development or results review, contact us to discuss how we can help.



Data Import

To help make the transition to SHOEBOX Audiometry smooth, we can format and import your historic patient and employee demographic data. You'll be up and running in no time.



Our Story

We are SHOEBBOX Audiometry—makers of the world’s first easy-to-use, automated, iPad audiometer. Why did we invent an audiometer on an iPad? Because we believe that everyone deserves access to hearing healthcare.

SHOEBBOX Audiometry ships in a small padded carry case and includes:

- ✓ Transducers
- ✓ Bone Oscillator and Switchbox (optional, Pro only)
- ✓ iPad stand

 www.shoebox.md

 Ottawa, ON, Canada

 info@shoebox.md

 1-877-349-9934 FREE (N.A.)
+1-613-728-6666 (Int'l)
+1-613-627-4798 (Fax)