


SHOEBOX QuickTest

A Simple Solution for New Customer Acquisition

Qualified leads for your hearing services clinic, at a fraction of the cost of traditional marketing.

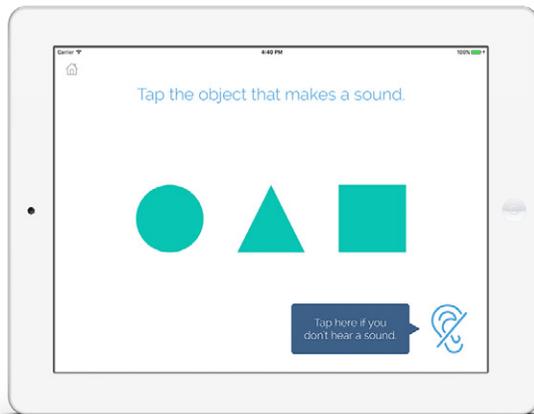

- ✓ Quick & simple screening test
- ✓ Flexible test options
- ✓ Custom branding
- ✓ Ability to manage a fleet of iPads through a web portal
- ✓ Easy and secure access to client test results
- ✓ Flexible subscription to minimize upfront costs

A challenge we often hear from our customers who own or operate hearing clinics is how best to attract qualified clients for audiology appointments and hearing aid fittings. SHOEBBOX® QuickTest is an easy-to-use and cost-effective new customer acquisition tool.

Use the system in attended or unattended kiosk mode at health fairs, in retail locations, waiting rooms, or almost any setting where your potential clients may gather. This quick, self-administered hearing test functions like playing a game. It can be configured to test up to four frequencies and the volume settings can be adjusted to help control for various testing environments. The user provides contact details—which can be customized to your specific needs—such as email address and phone number before being presented with their screening results. These results are automatically pushed to a cloud-based portal so you can quickly and easily identify your best candidates for scheduling follow-up appointments.

SHOEBBOX QuickTest is the ultimate companion to our Pro edition which can be used to perform a diagnostic hearing test after using QuickTest.



Easy-to-use

Simple, self-administered hearing screening test



Qualified Leads

Prioritize how you follow up based on user results



Cost Effective

Place in high traffic areas and watch the leads roll in

SHOEBOX QuickTest Features

Equipment Compatible with most iPads, contact us for more details

Recommended Transducers RadioEar DD450

Functionality

- Fast and fully automated screening test (1 minute)
- Easy-to-use game interface
- Up to 4 configurable frequencies
- Configurable upper primer and minimum tone output levels
- Collection of client information with a standard survey
- Surveys can be tailored through the optional customization service
- Customization of display with your own logo
- Centralized management of client screening results and survey information from all devices via the web portal

Sound Generation Range: 500Hz - 8000Hz
10 dB to 80 dB

Subscription Includes:

- ✓ Software and license to use on one iPad device
- ✓ Transducers
- ✓ Cloud-based web portal
- ✓ One hour training/setup session

Additional Options:

- ✓ Use your own iPad
- ✓ Survey and branding customization services
- ✓ Kiosk stand
- ✓ Remote iPad Management